

**HF1970 COMFY HYDRAULIC STRETCHER TROLLEY, 5 Functions**



**Standard Features:**

- > COMFY hydraulic stretcher trolley with adjustable back rest, knee rest, trendelenburg, reverse trendelenburg & adjustable height.
- > Back Rest adjustable by gas spring, knee rest adjustable by super smooth crank, other functions adjustable by hydraulic mechanism.
- > Size: 2082 × 860 × 608/923 mm.
- > Framework made of mild steel epoxy coated finish with X-ray permeable Top.
- > Collapsible side railings.
- > Equipped with hydraulic pump
- > Equipped with I.V. Pole, drainage hook & bumper protection.
- > 150mm dia noiseless castors with simultaneous braking system which locks / unlocks 2 castors with single pedal press.
- > International quality structure and functionality to satisfy the needs of the best hospital and Operation theaters

**HF1971** : Same as HF1970 but with two section platform, back rest is X-ray permeable top and wire mesh body section.


## HF1974 COMFY STRETCHER TROLLEY, 4 Functions



### Standard Features:

- > COMFY hydraulic stretcher trolley with adjustable back rest, trendelenburg, reverse trendelenburg & adjustable height.
- > Back Rest adjustable by gas spring, other functions adjustable by hydraulic mechanism.
- > Size: 1930 × 688 × 567/797 mm.
- > Equipped with hydraulic pump.
- > Equipped with I.V. Pole, drainage hook & bumper protection.
- > 150mm dia noiseless castors with simultaneous braking system which locks / unlocks 2 castors with single pedal press.
- > Two easy lifting guardrails (One on each side) made of PE engineering plastic, which are safe and reliable, and can be fixed upward and downward.
- > International quality structure and functionality to satisfy the needs of the best hospital and Operation theaters

## HF1975 HYDRAULIC STRETCHER TROLLEY, 2 Functions

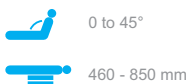

### Standard Features:

- > Hydraulic stretcher trolley with adjustable Back rest & adjustable height.
- > Back Rest adjustable by gas spring & Height adjustable by hydraulic mechanism.
- > Size: 2100 × 690 × 460/850 mm.
- > Framework made of mild steel epoxy coated finish with mesh mattress platform.
- > Collapsible side railings.
- > Equipped with bumper protection.
- > 125mm dia noiseless castors with simultaneous braking system which locks / unlocks 2 castors with single pedal press.
- > International quality structure and functionality to satisfy the needs of the best hospital and Operation theaters

### HF1977 STRETCHER TROLLEY, 2 Functions



 0 to 45°

 530 - 820 mm

#### Standard Features:

- > Stretcher trolley with adjustable Back rest & adjustable height.
- > Back Rest adjustable by gas spring & Height adjustable by crank mechanism.
- > Size: 1930 × 595 × 530/820 mm.
- > Framework made of mild steel epoxy coated finish with ABS Plastic Top.
- > Two easy lifting guardrails (One on each side) made of PE engineering plastic, which are safe and reliable, and can be fixed upward and downward.
- > 125mm dia noiseless castors with simultaneous braking system which locks / unlocks 2 castors with single pedal press.
- > Equipped with separate directional wheel (5th Wheel) for easy operation only by one person.
- > International quality structure and functionality to satisfy the needs of the best hospital and Operation theaters

### HF1980 STRETCHER TROLLEY

#### Standard Features:

- > Size: 1950L × 650W × 680 / 900 H mm.
- > Height adjustable from 600 - 950 H mm.
- > Frame work made of stainless steel.
- > S.S. Safety guard rails on each side.
- > Mounted on four castor (all with individual brakes).
- > Provision for IV rod.
- > Supplied completely with S.S. side railings and IV rod.



### HF1982 STRETCHER TROLLEY



#### Standard Features:

- > Size: 1950L × 650W × 780H mm.
- > Frame work made of stainless steel.
- > S.S. Safety guard rails on each side.
- > Basket for keeping blankets etc.
- > Mounted on four castor.
- > Provision for IV rod.
- > Supplied completely with S.S. side railings and IV rod.
- > With removable top

# □ MANUAL TRANSFER STRETCHER



## HF1981 PATIENT TRANSFER STRETCHER, 2 FUNCTION



### Standard Features:

- > Size: 1940 x 680 x 610-920mm
- > Height adjustable by crank mechanism (920 mm-610 mm)
- > Parallel Mattress surface movement by crank to transfer patient to Bed or table.
- > Horizontal movement range : 400mm
- > Aluminum – Alloy Guard rails, can be fixed up and down.
- > 150mm Central – locking castors, Equipped with separate directional wheel for easy maneuvering.


## HF1988 TRANSFER STRETCHER CUM MEDICAL TROLLEY, STAINLESS STEEL, Adjustable Height

### Standard Features:

- > Stainless Steel construction, for long life.
- > Size: 3500 × 650 × 750 / 950 mm.
- > Adjustable height from 750 - 950 mm by unobtrusive crank mechanism. Adjustable back section.
- > The two stretcher base can move up / down independently and separately, and can be brought to the same level for smooth slide.
- > Both base are locked securely before the bed slides.
- > All wheels with individual brakes.



## HF1989 COMFY TRANSFER STRETCHER CUM MEDICAL TROLLEY

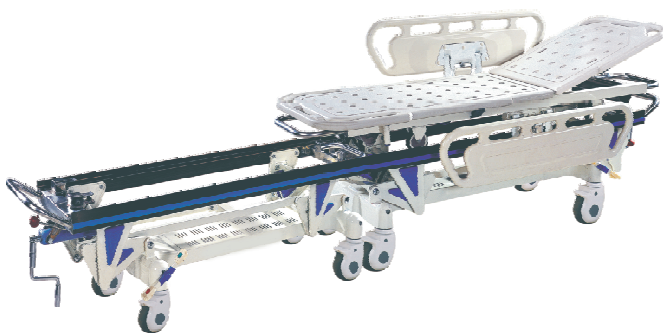

### Standard Features:

- > Made of high quality Aluminum alloy, PE coated, die pressed using automatic power press machines.
- > Size: 3650 × 640 × 640 / 970 mm. Adjustable back section.
- > Adjustable height from 640 - 970 mm.
- > Super strong aluminum alloy tracks, increase the durability, steadiness and reliability of the moving stretcher.
- > The two single-stretchers can move up / down independently, and can be conveniently brought to the same level for smooth slide.
- > Both base are locked securely and automatically before the bed slides.
- > Both single-stretchers are equipped with separate directional wheel (5th Wheel) for easy operation only by one person
- > 125mm dia noiseless castors with simultaneous braking system which locks / unlocks 2 castors with single pedal press.
- > International quality structure and functionality to satisfy the needs of the best hospital and Operation theaters.